DEL MEDICAL

Elevating, Four-Way Float Top Table

The Del Medical EV800, with an industry first 800lb patient weight limit, is built to demanding hospital quality standards and is ideally suited for any hospital or busy imaging center

The EV800's most notable attributes are its rack and pinion motor drives providing smooth, quiet table elevation and its precision roller bearing system that allows our floating table top to glide freely.

The smooth table top surface and soft rounded corners make the EV800 an attractive elevating table that is sure to compliment any imaging suite.

The EV800 is configurable with any Del Medical tube and/or wall stand for a complete imaging suite.

Product Features

- Smooth, quiet table operation with rack and pinion motor drive
- Integral collision sensors
- Heavy-duty electric table locks with foot treadle release
- All steel structure and welded base provides maximum strength and reliability
- Smooth, stain resistant fiber resin top





EV800

Elevating, Four-Way Float Top Table



Technical Specifications

Load capacity: 800lb (363 kg)
Base length: 53" (135cm)
Base width: 25" (64cm)

Table top movement: Longitudinal ± 21.25" (54cm)

Transverse ± 4.5 " (11.4cm)

Variable height: 22.25" (56.5cm) to 34.25" (87cm)

Bucky travel: ± 8.5 (22cm)

Some flat panel DR detectors may

have some affect on table

bucky travel.

Top length, standard: -optional 86.5" (220cm) 78" (198cm)

Top width, standard: -optional 31.9" (81cm) 35.9" (91cm)

Table weight: 563lb (256kg)

Aluminum equivalency: 0.7mm

Power Requirements: 120 VAC or 230-240 VAC

Certifications

- ETL listed


Product Options

- 17" x 17" (43cm x 43cm) grid cabinet or bucky
- Deluxe, heavy-duty manual cassette tray
- Size-sensing cassette tray with PBL interface
- Ion chamber (AEC)
- 78" (198cm) table top length
- 35.9" (91cm) table top width
- Carbon fiber composite tabletop, (0.4mm Al. eq.@ 100kVp). Available for standard width only.
- Variety of grid options
- · Lateral cassette holder
- Compression device
- · Patient hand grips







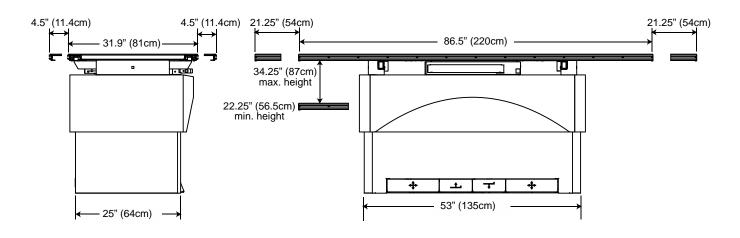

DEL MEDICAL

Toll Free: (800) 800-6006 Fax: (847) 288-7011 Email: info@delmedical.com www.delmedical.com Del Medical is a UMG Company. Del Medical is an ISO 13485 & ISO9001 Certified Facility. Specifications are subject to change without notice. Toll Free: (800) 261-9808 Fax: (914) 835-6111 Email: info@umgxray.com

www.umgxray.com 28 Calvert Street, Harrison, NY 10528,